Ab-DeliverIN™
Antibody Delivery Inside Living Cells

To find the ideal conditions, Ab-DeliverIN™ must be tested at ratios of **1 µL/µg**, **2 µL/µg**, and **2.5 µL/µg** (µL of Ab-DeliverIN / µg of Ab). For the antibody quantity, we suggest **0.4 µg** per well in 96-well, **1 µg** per well in 24-well and **5 µg** per well in 6-well.

1. Prepare 3 wells: seed the cells to be at 70% confluent on the day of transfection.



2. Prepare Antibody solution at 100 µg/mL in PBS

The presence of BSA as additives in antibody reagent can completely inhibit the antibody delivery. In case of BSA, remove it before proceeding.



3. Prepare 3 tubes of Ab-DeliverIN



96-well plate	24-well plate	6-well plate
0.4 µL/0.8 µL/1 µL	1 µL/2 µL/2.5 µL	5 µL/10 µL/12.5 µL


4. Add the antibody solution to each tube of Ab-DeliverIN



Incubate 10 min at RT

96-well plate	24-well plate	6-well plate
0.4 µg (4 µL of 100 µg/mL per well)	1 µg (10 µL of 100 µg/mL per well)	5 µg (50 µL of 100 µg/mL per well)

5. Add serum-free medium to each tube and distribute to cells

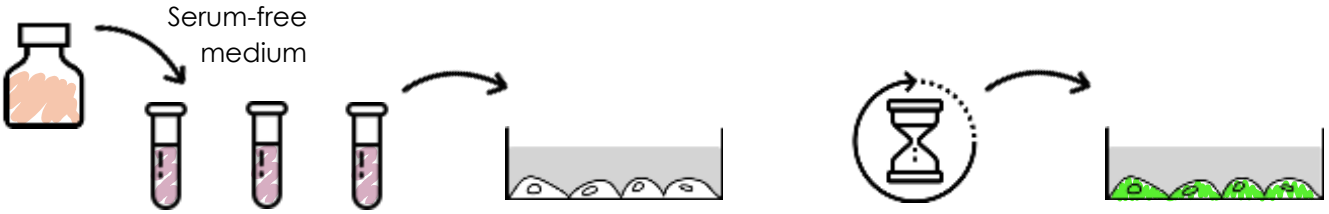

Up to 20 µL (96-well), 100 µL (24-well) or 200 µL (6-well).	Distribute each mix dropwise onto cells	Incubate 3 to 48h at 37°C	Choose the best ratio Ab:Ab-DeliverIN

Contact us regarding optimization procedures or for any technical question at: **tech@ozbiosciences.com** or visit our website: **www.ozbiosciences.com**.

OZ Biosciences SAS | Ph: +33 (0) 486 948 516 | contact@ozbiosciences.com | **www.ozbiosciences.com**

OZ Biosciences USA Inc. | Ph: + 1-858-246-7840 | contactUSA@ozbiosciences.com

Ab-DeliverIN™
Antibody Delivery Inside Living Cells

Before you begin

Ab-DeliverIN™ is a lipopolyamine formulation specifically designed to achieve intracellular delivery of biologically active antibodies (Ab).

Delivery kinetics depend on cell type. Perform a first evaluation at 3 to 4 h, then follow up to 24 or 48 h if needed.

Protocol | General considerations

• Allow reagents to reach room temperature and gently vortex them before forming complexes.
• Depending on Ab properties, the amount tested can be doubled.
• FITC-labeled IgG is provided as a positive control: use 1 or 2 µL Ab-DeliverIN per µg antibody for delivery assay.

⚠ Do not dilute Ab-DeliverIN before use.
⚠ We recommend preparing the complexes in a 0.5 mL tube to keep working volumes low.
⚠ Prepare Antibody dilution at 100 µg/mL in PBS. Do not use culture medium for this step. Ab concentration may vary depending on the application.

1. Cells preparation:
One day before transfection prepare the cells according to table 1 (recommended 70% confluence).

Cell culture format	Adherent cell nb. / well	Suspension cell nb. / well
96-well	0.05-0.15 x 1.10^5	0.5-1.0 x 1.10^5
24-well	0.5-1 x 1.10^5	1.5-5 x 1.10^5
6-well	2-5 x 1.10^5	5-20 x 1.10^5

Table 1: Cell number suggested per well.

2. Preparation of Antibody/Ab-DeliverIN complexes:

a. Antibody: Prepare antibody solution at 100 µg/mL in PBS and add the desired quantity into an empty tube.

b. Ab-DeliverIN: Add 0.8 to 10 µL of Ab-DeliverIN to the antibody solution (refer to Table 2).

c. Gently mix by pipetting up & down and incubate the mixture at room temperature for **10 to 15 min**.

⚠ Do not vortex or centrifuge.

3. Transfection:
a. Add 20 to 200 µL serum-free medium to the antibody/Ab-DeliverIN mix (refer to table 2).

b. Disperse immediately drop-by-drop the complexes onto cells in their regular culture medium.

c. Gently rock the plate and cultivate the cells at **37°C** in a CO_2 incubator under standard conditions until evaluation of delivery efficiency (3 to 48h).

Cell culture format	Antibody quantity	Ab-DeliverIN volume	Dilution volume	Transfection volume
96-well	0.4 µg	0.4 µL	20 µL	200 µL
24-well	1 µg	2 µL	100 µL	500 µL
6-well	5 µg	10 µL	200 µL	2 mL

Table 2: Recommended Antibody amounts, Ab-DeliverIN volumes and transfection conditions. For other formats, adjust conditions accordingly.

Protocol Optimization

Contact us regarding protocols for optimization procedures, or any technical question:
✉ **- tech@ozbiosciences.com.**

Use, handling and storage

For Research Use Only. Not for use in humans. Not for use in diagnostic or therapeutic purposes.

Shipping conditions: Room Temperature
Storage conditions Ab-DeliverIN: +4°C
Shelf life: 1 year from the date of purchase

Kit contents

AI20100 : 100 µL of Ab-DeliverIN reagent
AI20250 : 250 µL of Ab-DeliverIN reagent
AI20500 : 500 µL of Ab-DeliverIN reagent
AI21000 : 1 mL of Ab-DeliverIN reagent

Ab-DeliverIN is provided with 100 µL of FITC-labeled IgG positive control.

Certificate of analysis on demand.

Related Products

Ref	Description
#PI11000	Pro-DeliverIN protein delivery reagent
#IGG0100	IgG-FITC (100 µg/mL) in sodium azide pH 7.4

Purchaser Notification | Conditions of Sale

This product is sold in accordance with our general conditions of sale that you can find on our website: https://ozbiosciences.com/content/3-terms-and-conditions.

OZ Biosciences SAS | Ph: +33 (0) 486 948 516 | contact@ozbiosciences.com | **www.ozbiosciences.com**

OZ Biosciences USA Inc. | Ph: + 1-858-246-7840 | contactUSA@ozbiosciences.com