

## pVectOZ-LacZ ( $\beta$ -Galactosidase Expression Vector)

### Description

The pVectOZ-LacZ vector has been created to produce the highest levels of  $\beta$ -Galactosidase expression in a broad range of mammalian cells and tissues. It contains a proprietary modified human cytomegalovirus (CMV) promoter followed by a specific intron, enhancer and a terminator. The expression vector is engineered in an optimized plasmid backbone to achieve the highest levels of transgene expression in mammalian cells and high copy number production in *Escherichia coli*.

### Kit contents

Ref. #PL00030  
25  $\mu$ g pVectOZ-LacZ (encoding for  $\beta$ -Galactosidase) plasmid in 25  $\mu$ l sterile TE buffer.

### Storage

Store at -20°C.

### Selection Marker

**Kanamycin** is the selection gene included for producing the plasmid in *Escherichia coli*.

### Applications

pVectOZ-LacZ ( $\beta$ -Galactosidase) vector is suitable for all transfection applications (in vitro & in vivo).

**Presentation.** The transgene expression level depends mainly on the promoter, enhancer, terminator and plasmid backbone. The pVectOZ-LacZ expression cassette was designed to express very high levels of transgene product in many mammalian cells and tissues. This vector has been modified to eliminate sequences affecting transgene expression levels while optimizing those critical for high levels of expression. The final expression cassette accommodates the high levels of transgene expression in mammalian cells as well as high yield of plasmid production in *Escherichia coli*. The resulting plasmid is the ideal vector to reach high levels of expression in vitro and in vivo.

**Use.** For high levels of transgene expression in mammalian cells and tissues. For optimal results, this vector can be used with all OZ Biosciences transfection reagents to transfect a wide variety of mammalian cells and tissues.

### $\beta$ -Galactosidase detection

LacZ is one of the most frequently used reporter gene in transfection experiments because of the gene product specific properties. Indeed, the LacZ encoded protein,  $\beta$ -galactosidase, is very stable, resistant to proteolytic degradation and easily tested. The levels of active  $\beta$ -galactosidase expression can be quickly measured by colorimetric assays:

- CPRG (Chlorophenol red- $\beta$ -D-galactopyranoside) assay kit (catalog # GC10002)
- ONPG (o-nitrophenyl- $\beta$ -D-galactopyranoside) assay kit (catalog # GO10001)

The levels of active  $\beta$ -galactosidase expression can also be detected by histochemistry:

- X-gal staining kit (catalog # GX10003). It allows to visualize  $\beta$ -gal expression, through hydrolysis of X-Gal which yields a blue precipitate. Thus, transfected cells appear blue and can be visualized by microscopy.

### References

1. Felgner, J. *et al.* (1994). *J. Biol. Chem.* **269**: 2550.
2. Gussoni, E. *et al.* (1996). *Nature Biotech.* **14**: 1012.
3. Liu, M. *et al.* (2003). *Methods Mol Biol.* **235**: 289-96.
4. Nazarenko, DA. *et al.* (2001). *Biotechniques*, **4**: 776.
5. Tanahashi, H. *et al.* (2000). *Anal Biochem.* **279**: 122.

### Contact Us


Please, feel free to contact us for all complementary information and remember to visit our website ([www.ozbiosciences.com](http://www.ozbiosciences.com)) to stay informed on the latest breakthrough technologies and updates on our complete product list.

Technical questions: [tech@ozbiosciences.com](mailto:tech@ozbiosciences.com)  
Order: [order@ozbiosciences.com](mailto:order@ozbiosciences.com)  
General Information: [contact@ozbiosciences.com](mailto:contact@ozbiosciences.com)

**OZ Biosciences SAS**  
163 avenue de Luminy  
Case 922, zone entreprise  
1328 Marseille cedex 09 - FRANCE  
Ph: +33 (0) 486 948 516  
Fax: +33 (0) 463 740 015 [c](http://www.ozbiosciences.com)  
[contact@ozbiosciences.com](mailto:contact@ozbiosciences.com)  
[order@ozbiosciences.com](mailto:order@ozbiosciences.com)

**OZ Biosciences INC**  
7975 Dunbrook Road,  
Suite B  
San Diego CA 92126- USA  
Ph : + 1-858-246-7840  
Fax : + 1-855-631-0626  
[contactUSA@ozbiosciences.com](mailto:contactUSA@ozbiosciences.com)  
[orderUSA@ozbiosciences.com](mailto:orderUSA@ozbiosciences.com)

## Plasmid Map

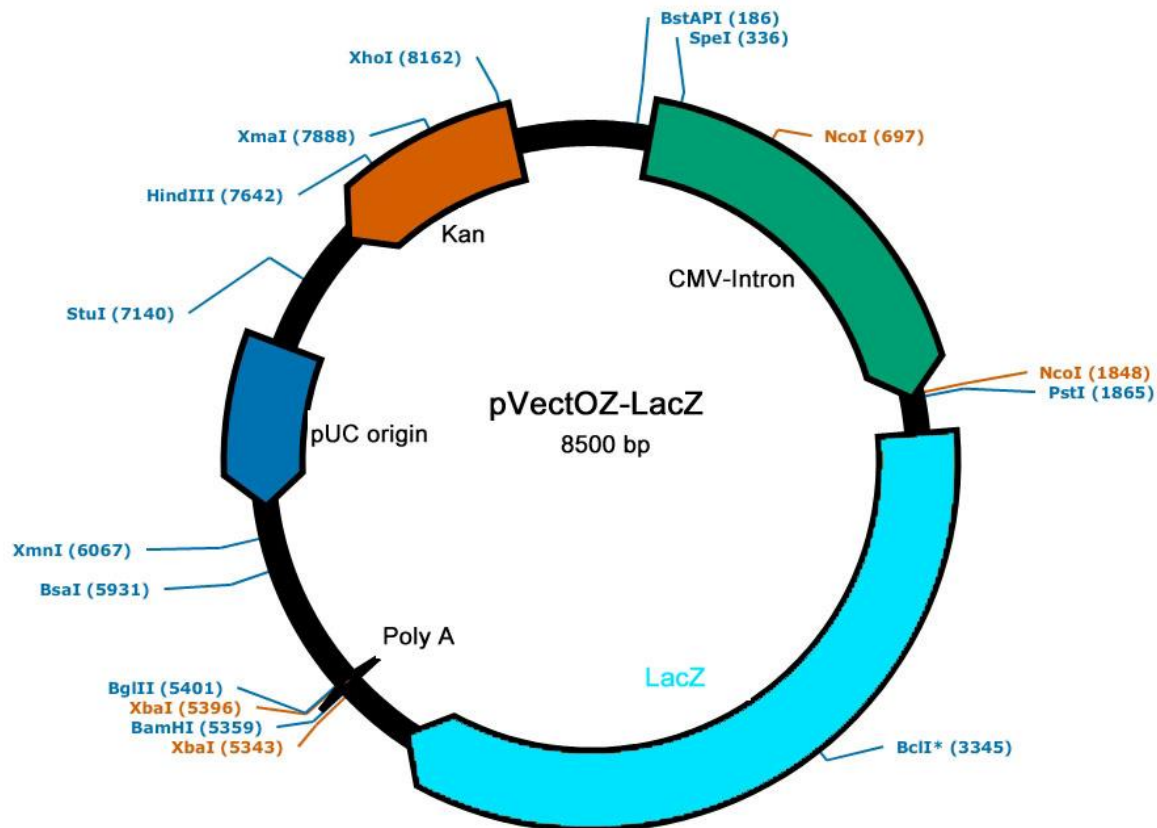

## Related Products

Description	Reference	Description	Reference
<b>Magnetofection Technology</b>		<b>Gene &amp; Protein Tools</b>	
Mega Magnetic Plate	MF14000	Bradford – Protein Assay Kit	BA00100
Super Magnetic Plate	MF10000	GeneBlaster selection kit	GB20010
Magnetic Plate 96-magnets	MF10096	β-Galactosidase (ONPG) assay kit	GO10001
PolyMag 1mL (for all nucleic acids)	PN31000	β-Galactosidase (CPRG) assay kit	GC10002
PolyMag Neo 1mL (for all nucleic acids)	PG61000	X-Gal Staining Kit	GX10003
CombiMag 1mL (boost transfection reagent)	CM21000	<b>Plasmids</b>	
SilenceMag 500μL (for siRNA applications)	SM10500	pVectOZ-CAT 25μg	PL00010
NeuroMag 1mL (for neuron transfection)	NM51000	pVectOZ-GFP 25μg	PL00020
		pVectOZ-Luc 25μg	PL00040
		pVectOZ-SEAP 25μg	PL00050
<b>Lipofection (lipid-based reagents)</b>		pVectOZ-CAT 100μg	PL00110
DreamFect Gold Transfection reagent 1mL	DG81000	pVectOZ-GFP 100μg	PL00120
DreamFect Transfection reagent 1mL	DF41000	pVectOZ-LacZ 100μg	PL00130
Lullaby siRNA Transfection reagent 1mL	LL71000	pVectOZ-Luc 100μg	PL00140
VeroFect Transfection Reagent 1mL	VF61000	pVectOZ-SEAP100μg	PL00150
FlyFectin Transfection Reagent 1mL	FF51000		

Rev 10/20 FB