

pVectOZ-Luc (Luciferase Expression Vector)

Description

The pVectOZ-Luc vector has been created to produce the highest levels of Luciferase expression in a broad range of mammalian cells and tissues. It contains a proprietary modified human cytomegalovirus (CMV) promoter followed by a specific intron, enhancer and a terminator. The expression vector is engineered in an optimized plasmid backbone to achieve the highest levels of transgene expression in mammalian cells and high copy number production in *Escherichia coli*.

Kit contents

Ref. #PL00140
100µg pVectOZ-Luc (encoding for Secreted Alkaline Phosphatase) plasmid in 100µl sterile TE buffer.

Storage

Store at -20°C.

Selection Marker

Kanamycin is the selection gene included for producing the plasmid in *Escherichia coli*.

Applications

pVectOZ-Luc (Luciferase) vector is suitable for all transfection applications (*in vitro* & *in vivo*).

Presentation. The transgene expression level depends mainly on the promoter, enhancer, terminator and plasmid backbone. The pVectOZ-Luc expression cassette was designed to express very high levels of transgene product in many mammalian cells and tissues. This vector has been modified to eliminate sequences affecting transgene expression levels while optimizing those critical for high levels of expression. The final expression cassette accommodates the high levels of transgene expression in mammalian cells as well as high yield of plasmid production in *Escherichia coli*. The resulting plasmid is the ideal vector to reach high levels of expression *in vitro* and *in vivo*.

Use. For high levels of transgene expression in mammalian cells and tissues. For optimal results, this vector can be used with all OZ Biosciences transfection reagents to transfect a wide variety of mammalian cells and tissues.

Luciferase detection

For transfections performed with pVectOZ-Luc plasmids, the detection can be realized by using standard Luciferase assay kits.

References

1. Bronstein, I. *et al.* (1994). *Anal Biochem.* **219**: 169.
2. Welsh, S. *et al.* (1997). *Curr Opin Biotechnol.* **8**: 617.
3. Doshida, M. *et al.* (2006). *J Biol Chem.* **281**: 24270.
4. Deleuze, V. *et al.* (2007). *Mol Cell Biol.* **27**: 2687-97.
5. Mykhaylyk, O. *et al.* (2008). *COMT.* **10**: 493.
6. Luker, KE. *et al.* (2008). *Antiviral Res.* **78**: 179.
7. Burdyga, G. *et al.* (2008). *J Neurosc.* **28**: 11583.
8. Arnoux, V. *et al.* (2008). *Mol Cell Biol.* **19**: 4738.

Contact Us


Please, feel free to contact us for all complementary information and remember to visit our website (www.ozbiosciences.com) to stay informed on the latest breakthrough technologies and updated on our complete product list.

Technical questions: tech@ozbiosciences.com
Order: order@ozbiosciences.com
General Information: contact@ozbiosciences.com

OZ Biosciences SAS
163 avenue de Luminy
Case 922, zone entreprise
13288 Marseille cedex 09 - FRANCE
Ph: +33 (0) 486 948 516
Fax: +33 (0) 463 740 015
contact@ozbiosciences.com
order@ozbiosciences.com

OZ Biosciences INC
7975 Dunbrook Road,
Suite B
San Diego CA 92126- USA
Ph : + 1-858-246-7840
Fax : + 1-855-631-0626
contactUSA@ozbiosciences.com
orderUSA@ozbiosciences.com

Plasmid Map

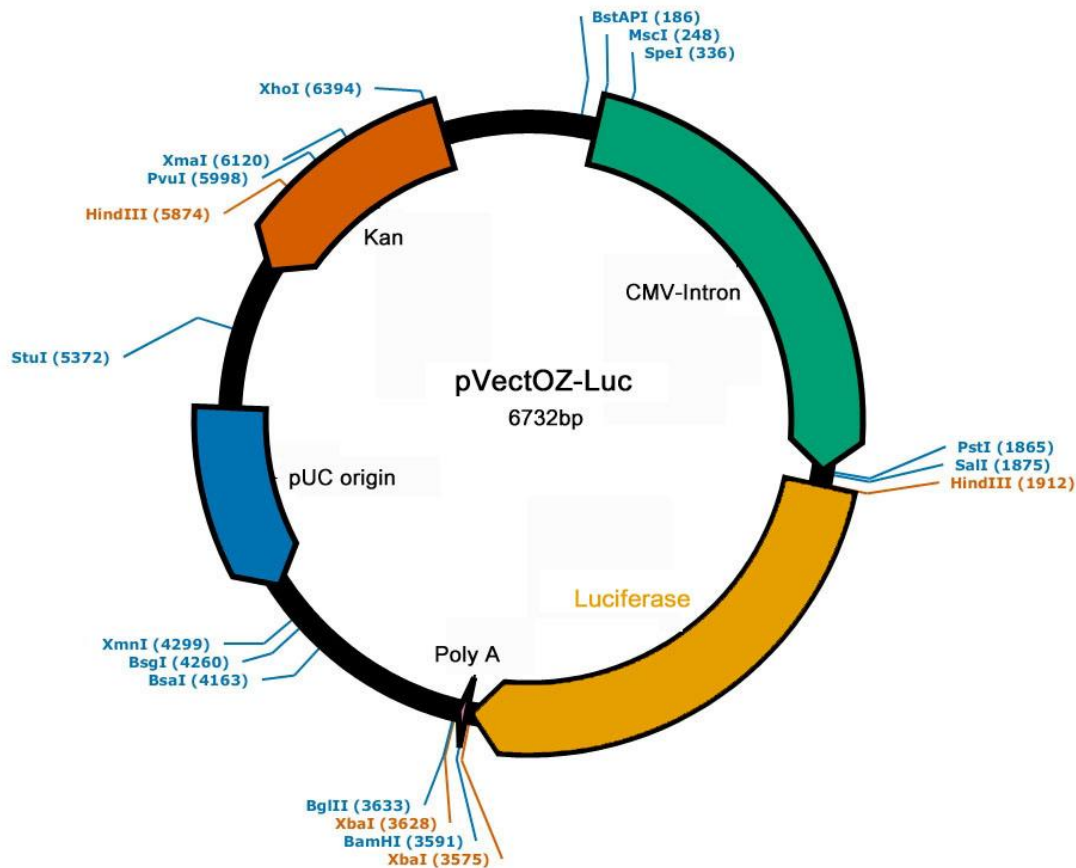

Related Products

Description	Reference	Description	Reference		
Magnetofection Technology		Gene & Protein Tools			
Mega Magnetic Plate	MF14000	Bradford – Protein Assay Kit	BA00100		
Super Magnetic Plate	MF10000	GeneBlaster selection kit	GB20010		
Magnetic Plate 96-magnets	MF10096	β-Galactosidase (ONPG) assay kit	GO10001		
PolyMag 1mL (for all nucleic acids)	PN31000	β-Galactosidase (CPRG) assay kit	GC10002		
PolyMag Neo 1mL (for all nucleic acids)	PG61000	X-Gal Staining Kit	GX10003		
CombiMag 1mL (boost transfection reagent)	CM21000	Plasmids			
SilenceMag 500μL (for siRNA applications)	SM10500	pVectOZ-CAT 25μg	PL00010		
NeuroMag 1mL (for neuron transfection)	NM51000	pVectOZ-GFP 25μg	PL00020		
Lipofection (lipid-based reagents)		pVectOZ-LacZ 25μg	PL00030		
		DreamFect Gold Transfection reagent 1mL	DG81000	pVectOZ-Luc 25μg	PL00040
		DreamFect Transfection reagent 1mL	DF41000	pVectOZ-SEAP 25μg	PL00050
		Lullaby siRNA Transfection reagent 1mL	LL71000	pVectOZ-CAT 100μg	PL00110
		VeroFect Transfection Reagent 1mL	VF61000	pVectOZ-GFP 100μg	PL00120
		FlyFectin Transfection Reagent 1mL	FF51000	pVectOZ-LacZ 100μg	PL00130
				pVectOZ-SEAP100μg	PL00150

Rev 10/20 FB